



www.flatmedical.com

US

(+1) 878-313-4193
sales@flatmed.org
1423 Broadway Ste 214,
Oakland, CA 94612, USA

EU

(+31) 970-1020-474
europe@flatmed.org
Chromiumweg 139, 3812 NM,
Amersfoort, Netherlands


TW

(+886) 2-2567-2959
info@flatmedical.com
12F-2, No. 206, Sec. 2, Nanjing E. Rd.,
Zhongshan Dist., Taipei City 104496, Taiwan



EpiFaith[®] Syringe


A syringe designed for epidural space localization with a confidence signal

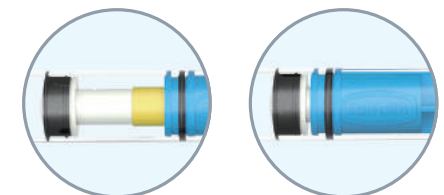

EpiFaith[®] Syringe

Hassle-free and reliable design that promotes safer epidural injections.

EpiFaith[®] Syringe, a compact and sturdy device, assists physicians in locating the epidural space. Intuitively equipped with the visual confidence signal upon entry into the epidural space. The reassuring nature of the EpiFaith[®] Syringe has made the process safer and more convenient for both junior and senior physicians.

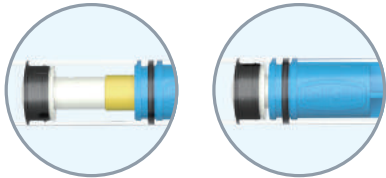

LOR Mechanism Explained:
 The pressure is proportionally increased as the distance between the plunger and the rubber piston gets closer. As the needle tip enters the epidural space, a rapid decrease in pressure also known as 'loss of resistance', causes the swift movement of the colour indicator ring to appear, signalling the endpoint for the advancement of the needle.



CLEAR SIGNAL

Provides real-time objective and quantitative pressure information, triggering a visual signal upon entry into the epidural space.



INTUITIVE CONTROL

Easy push/pull piston control allows physicians to monitor pressure while using both hands to advance the needle.




ADAPTABILITY

Compatible with other mainstream products used in epidural procedure and current practices.



Order Information

Item name	Type	Ref.#	Qty/Case
EpiFaith Syringe	Luer slip (US)	FM-01SLR	50
	NRFit (US)	FM-02SNR	
	Luer slip (EU)	EF-01SLR	
	NRFit (EU)	EF-02SNR	

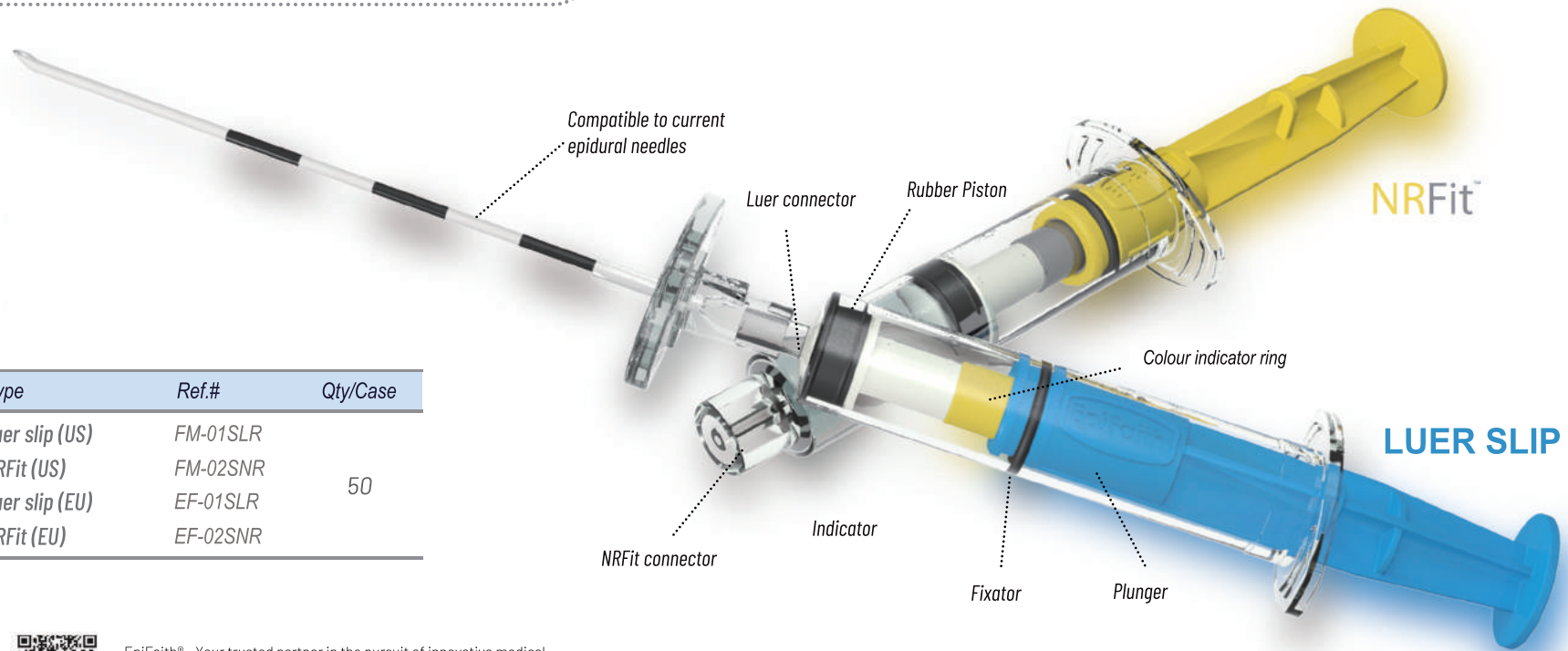

EpiFaith[®] - Your trusted partner in the pursuit of innovative medical solutions. Reach out today for more information or to place an order!
www.flatmedical.com